



Electric Music Toothbrush Operation Manual

Please Read It Carefully Before Using

www.osmiler.com

Shenzhen Osmiler Technology Co., Ltd.

Hello

First of all, thank you for choosing our sonic electric toothbrush. Let us have the opportunity to provide you with our high-quality service.

In order to bring you the enjoyment of high-quality service. After purchasing this product, please read this instruction carefully, and put the warranty card away prudently.

Product List


Brush head + transparent toothbrush cover	×2
The toothbrush body	×1
Type-C data cable	×1
Operating manual/warranty card	×1

Please keep the transparent toothbrush cover to protect the bristles on your toothbrush while traveling.


Matters needing attention:

Do not use the charger in the bath ! Do not change the battery at any time. If you need to replace the battery, please contact the manufacturer.

Product Figure

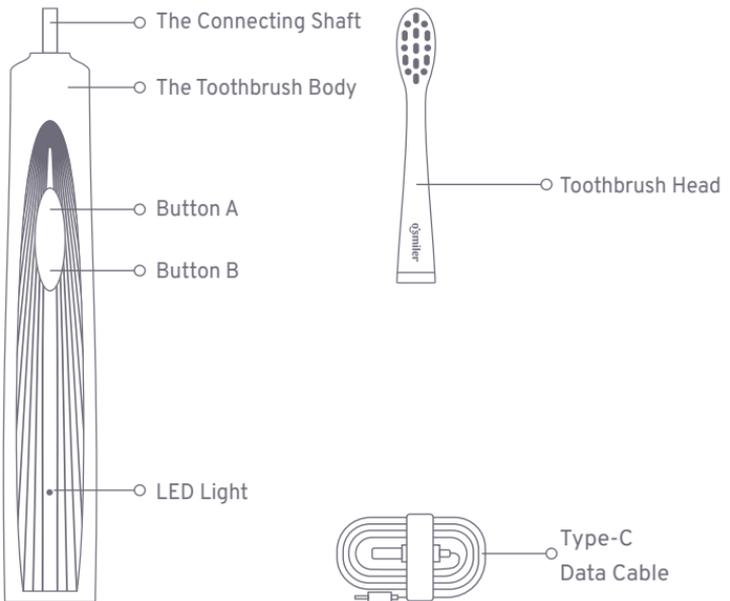

Function Graphic

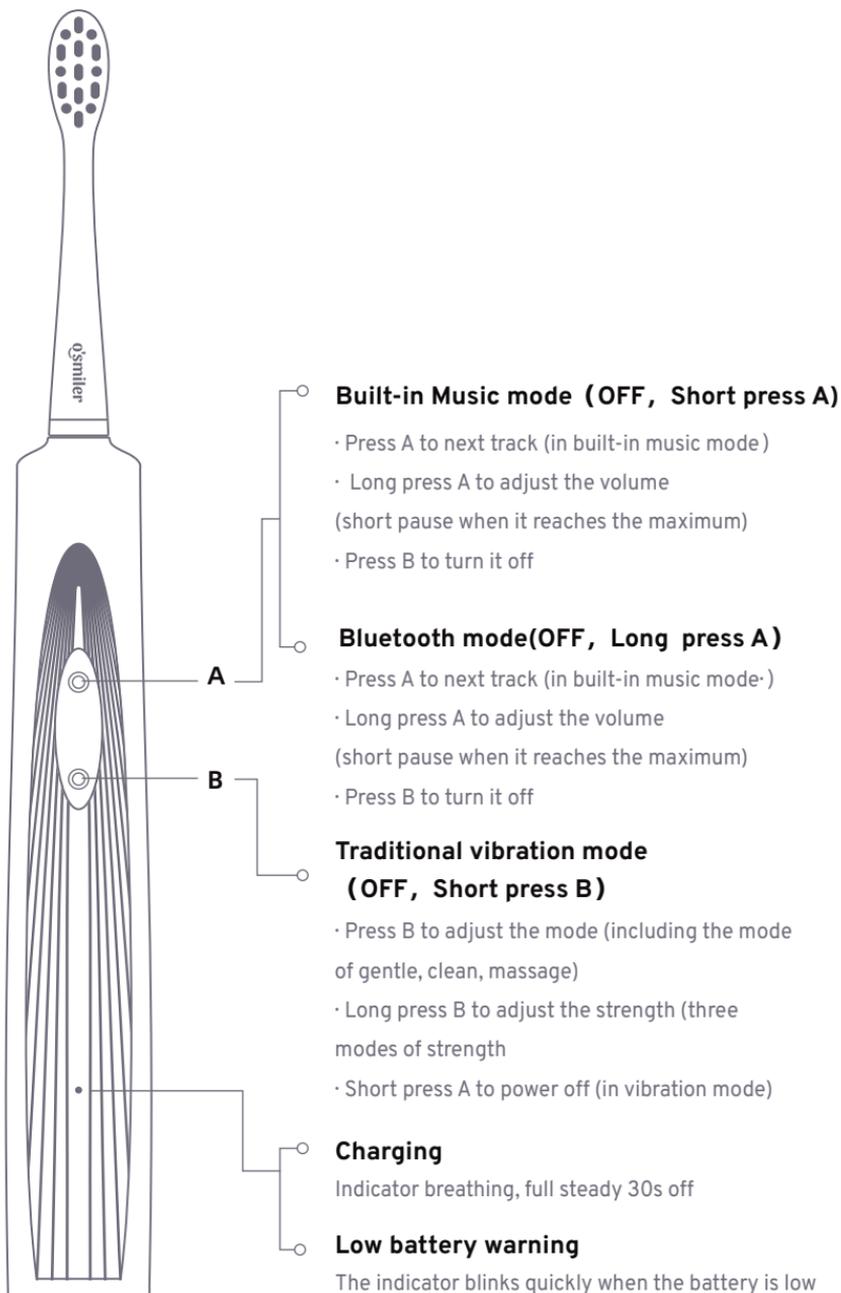

Assemble Electric Music Toothbrush

Insert the brush head downward into the connecting shaft of the body. There will be a gap (about 0.3-1.2mm) between the body and the brush head after insertion. This is the buffer distance reserved for the vibrating motor..





Charge

Connect the handle to the charging line and plug in the Type-c charger for charging. When charging, the host machine cannot be started. If charging is in operation, the host machine will be shut down;


- During charging, the indicator light is white for breathing
After being fully charged, the white light is steady on for 30 seconds and then off
- After being fully charged, the white light is steady on for 30 seconds and then off
- Charging time ≥ 4 hours

Upload Music


- 1 Connect the USB port of the type-C charging port at the bottom of the toothbrush with the type-C data cable attached in the package, and the indicator light is on
- 2 Once connected to the computer, the removable disk option (15.8 MB) will appear in the computer folder

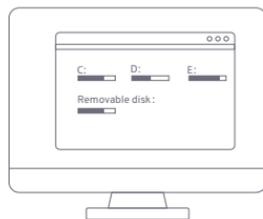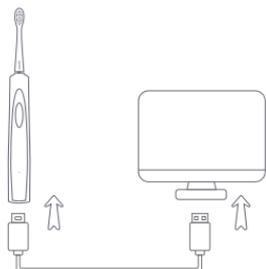

- 3 Copy selected mp3 format music to removable disk (keep music file size within 3-5mb)

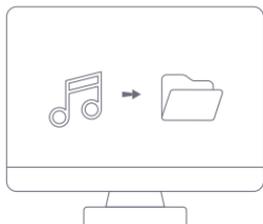


- 4 Roll out the removable disk on the computer, unplug the data line and use the product

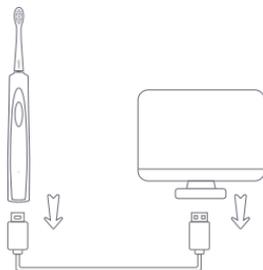

Built-in Music Mode


- 1 Apply appropriate amount of toothpaste to the bristles and soak the bristles

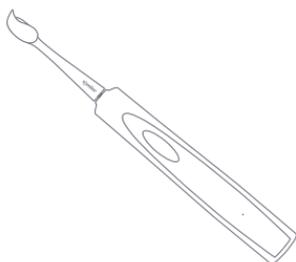


- 2 Put the brush head in the mouth, the bristles are close to the teeth, the toothbrush starts to work.




- 3 In the off state, Press A to enter the built-in music mode

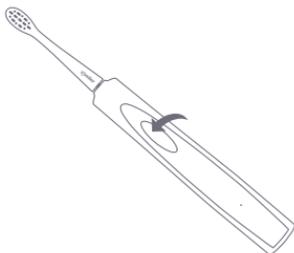


- 4 In built-in music mode, press A to switch songs

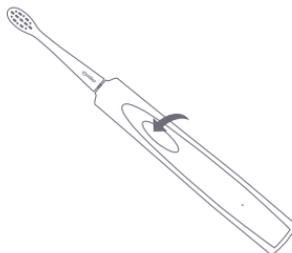

- 5 In built-in music mode, long press A to adjust to appropriate volume

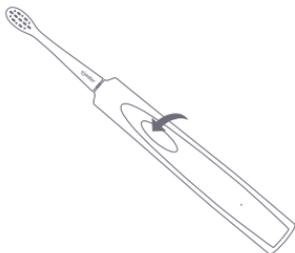


- 6 After brushing, press B to turn off and rinse with water Toothbrush body and brush head

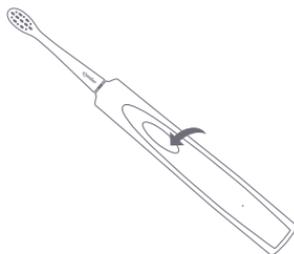

Bluetooth Mode


- 1 In the off state, long press A to enter the Bluetooth connection state. The toothbrush vibrates with Bluetooth connection prompt, and the indicator light flashes

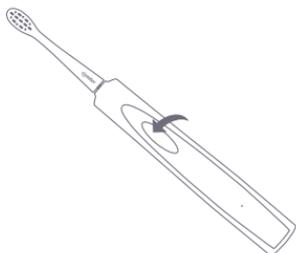


- 2 Open the Bluetooth interface of the mobile phone, select **Osmiler** as the Bluetooth name to connect, and the LED is always on



- 3 Open your favorite music player and select music playback.
- Short press A to switch music
 - Long press A to adjust the volume


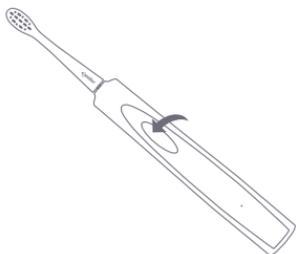


- 4 After brushing, press B to turn off and rinse with water Toothbrush body and brush head

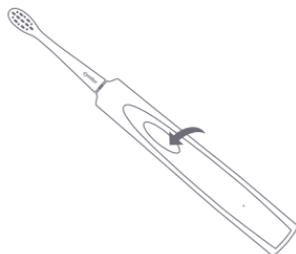


Traditional Vibration Mode

- 1 In the off state, press B to enter traditional vibration mode
- 2 In vibration mode, press B to select different Vibration mode (divided into three modes: gentle, clean and massage)

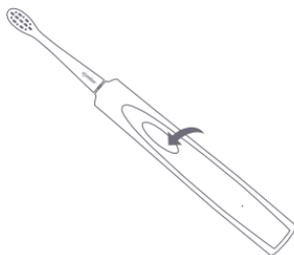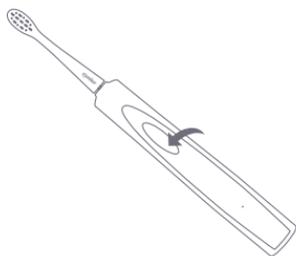

- 3 Under different vibration modes, long press B to adjust three different vibration forces
- 4 After brushing, briefly press A to turn off the machine.

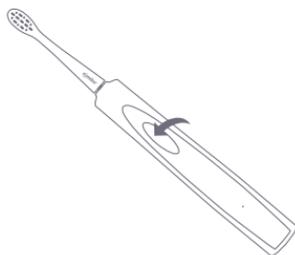

Tips for Brushing Teeth

- 1 For best teeth cleaning, change the brushing area every 30 seconds


Note:

To ensure that you brush your teeth evenly, please refer to the following four area


Outside of upper teeth




Inside of upper teeth

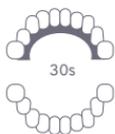


Outside of lower teeth

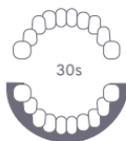

Inside of lower teeth

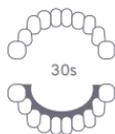


- 2 cleaning at 45° Angle, effect is better

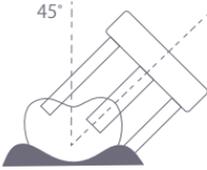


- 3 Do not press hard when brushing your teeth. Keeping the bristles and teeth intact have certain contact is enough




- 4 Remove and replace the brush head by pulling it up in the direction of the arrow (Do not twist the brush head)

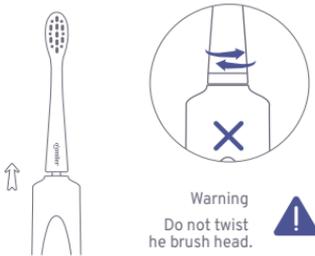


- 5 When you want to stop brushing and take your brush head out, Please turn off your toothbrush to avoid toothpaste splashing

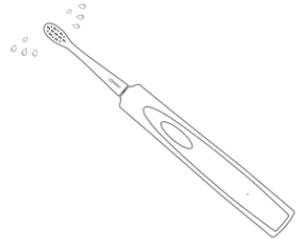

- 6 Keep the toothbrush in a dry and ventilated environment
Keep away from direct sunlight and fire and other heat-generating objects, otherwise it may lead to battery fluid leakage and fire and explosion, etc.

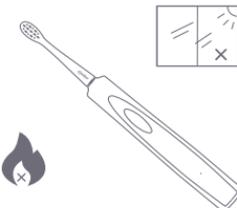


- 7 Rinse brush body and brush head directly with water. After rinsing, use soft paper towels or soft towels to wipe it clean

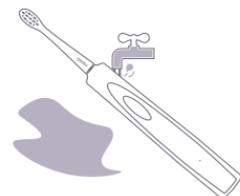

Product Parameters

Product model	S0224
Working current	500mA±100mA
Handle voltage rating	3.7V
Motor vibration frequency	250~330 Hz
Size	27.3 (L) X28.4 (W) X243.5mm(H)
Working power	3.0W
Waterproof level	IPX7
Weight	132.0g (with toothbrush head)
Battery	2000mAh
The toothbrush head	Food grade Dupont bristles with protective caps
Charging current	550±50mA
Certification	FCC, RoHS2.0, EMC, FDA

Weekly Maintenance

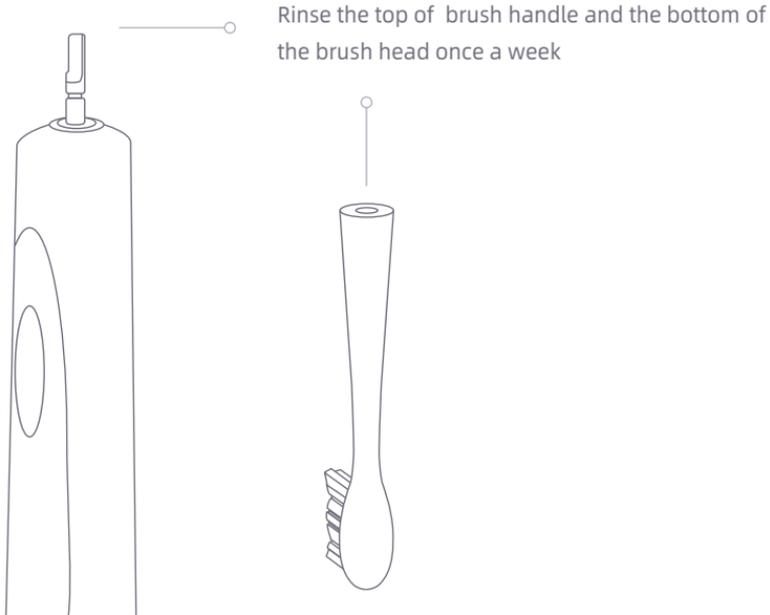

Security Service

Warning: Read all instructions before use

1. Check the product carefully before each use. If there is any damage to the product (brush head, handle, etc.), please stop using it.
2. This product does not contain parts that users can maintain by themselves. Do not disassemble, repair or modify the product. Only certified Professional and technical personnel are qualified to solve any problem. If the product no longer works properly or requires repair, Please read "Warranty And support " for reference
3. Please place and store products away from children or the disabled.
4. When not in use for a long time, please be sure to store the product in a cool and dry place, avoid direct sunlight.
5. In order to avoid damaging the product, do not put the brush head (including the protective cover) and handle into the dishwasher for cleaning.
6. Replace the Osmiler brush head every 3 months to work best. Replace the brush head with Osmiler brush head only.
7. Do not use toothbrush while charging.
8. Do not use the charging cable with damaged power cord.
9. Do not use in warm bath or hot shower. Although the product is IPX7 waterproof, it can not be effectively protected from hot water or steam.
10. This product is only used to clean teeth, gums and tongue. Use this product only for its intended use as described in this user manual or recommended by your dental professional.
11. This product is not suitable for the use of people (including children) with reduced physical, sensory or mental ability or lack of experience and knowledge. For their safety, please use it under the supervision of their responsible person or under the guidance of the use of electronic devices.
12. Children under 7 years old are not recommended to use this product.
13. This product is a personal care device and is not suitable for multiple patients in dental clinics or institutions.
14. When charging the toothbrush body, please use a UL listed computer, USB hub or USB duplex outlet.
15. To connect to a power source that is not in your country, use a correctly configured power outlet adapter.
16. If you have had oral or gum surgery in the past 2 months, please consult your dentist before using this product.
17. Contact your dentist if excessive bleeding occurs after use or if bleeding continues after using it 1 week.
18. Stop using this product and contact your doctor/dentist if you feel discomfort or pain.
19. This product complies with safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, please contact your doctor or device manufacturer before use.
20. If you have any other medical concerns, consult your doctor.
21. This toothbrush can be safely used with braces (the brush head wears out more quickly when used on braces) and restorations (fillers, crowns, veneers).

22. Itchy gums, numbness and even slight bleeding are normal when first used. Your gums will adjust to the sonic vibrations in about two days, and the discomfort will last for a week. If discomfort persists, consult your dentist.

23. If the product will not be used for a short period of time, please be sure to charge it every 6 months before use.

FCC Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

RF warning statement:

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Warranty and Support

QC PASS

Osmiler guarantees its products for one year from the date of purchase. Defects due to material and workmanship defects will be repaired or replaced at Osmiler's expense (including Osmiler retailers and dealers) if a convincing proof of purchase is provided during the warranty period.

Osmiler recommends using genuine Osmiler brush heads for optimal performance. Using other brush heads may damage your brush and limit your warranty service.

Contact your local retailer and dealer or email service@osmiler.com for warranty service.

Warranty and Support

Contents not covered by warranty:

- brush head.
- Damage caused by using unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, negligence, alteration or unauthorized repair.
- Normal wear and tear, including debris, scratches, scrapes, discoloration or fading.
- Proof of purchase has been altered, defaced or removed.

All implied warranties, including warranties of merchantability and fitness for a particular purpose, are limited to the duration of this warranty, nor any of its affiliates or subsidiaries shall be liable for any special, incidental, consequential or other damages based on breach of warranty, breach of contract, negligence, tort or any other legal theory. Such damages include, but are not limited to, lost savings or income, lost profits; loss of use; claims from third parties, including but not limited to dentists and dental hygienists; the cost of any alternative equipment or services. If the product needs to be returned to the retailer or dealer for repair or replacement, the following items must be provided for the repair and return of the product:

- 1.For warranty returns, include some form of proof of purchase, such as a mechanical replica or a duplicate sales receipt. Make sure the date of purchase and product model number are clearly legible. If the original has been sent, it cannot be returned;
- 2.Send the whole product;
- 3.Pack the product properly to prevent damage in transit. If possible, use original packing materials.
4. Attach a description of what happened to the product. Include the name and address of the location where the product is returned, printed or legibly printed, and a phone number (required for shipment). Please fill in the following customer information form:

Product model:

Name:

Purchase Date:

Contact information:

Email:

Address::

Problem Description: _____
